



i THE COMPRESSIVE CRYOTHERAPY

This recognized therapy combining ice water cryotherapy and air compression reduces the pain by two, accelerates healing and offers a faster recovery of physical activity.

The Moove sport will optimize your return on the field, it will be an incountournable help to your sporting deadlines and your ally for the recovery.

The sport Moove concerns the following indications:

Trauma management, muscle preparation and recovery, re-training, heavy legs.

Icy water passes through the splint, in combination with a continuous or alternating air pressure, so that the skin temperature goes from 30 ° to 6 ° C on average, which will relieve the treated area and speed up the recovery. limiting drug intake and promoting range of motion.



The Moove sport with its carrying bag and his knee splint

★ EFFECTS

- REDUCES** the pain
- DECREASES** the inflammatory process
- ABSORBS** edemas
- ACCELERATES** tissus scarring
- ACCELERATES** the recovery of activity

i THE CONCEPT OF EASYCRYO


Performance : The MOOVE model offers maintained pressure with the objective of treating deep tissue and alternating pressure to promote recovery and allow drainage of inflammation.

Efficiency : The intense cold will immediately relieve and block the inflammatory process. This technology is used in muscle preparation to increase endurance and recovery to drain toxins. 10 to 30 minutes of care be enough to relieve, revive exchanges, drain and relax the muscular area.

Simplicity : The ice water tank can also be cooled at freezer to extend the duration of use. Our devices operate on mains power supply and on battery too (5 hours) The set is easily portable. The device can be supplied with a carrying bag.



Moove sport

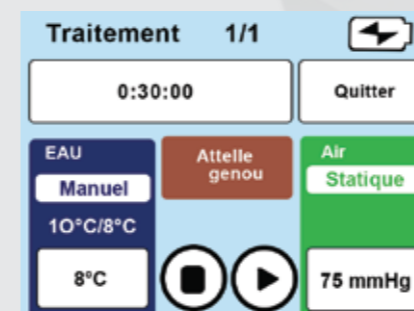

+ ADVANTAGES

- Programmable and modifiable system
- Automatic purging of splints for a better hygiene
- Machine washable splints protective covers
- Chilled water tank: with ice cubes or at the freezer (3h cold guarantee).



Water tank

THE SPLINTS



SHOULDER



THIGH



ANKLE

UPPER LIMB

- Shoulder right and left
- Elbow / wrist

LOWER LIMB

- Knee (L. 38 cm)
- Thigh (L. 50 cm)
- Ankle
- Hip right and left

TECHNICAL DETAILS

	MOOVE
Number of exits	1
min / max compression	0 mmHg / 75 mmHg
Compression	dynamic/constant
Water temperature	1° à 8°C recommended
Program	1 (alterable)
Cycle of Cold and Rest	Fonction repeat
CE med	LNE/GMED 0459


MOOVE SPORT

THE COMPRESSIVE CRYOTHERAPY

“

A professional tool for a personal use

”

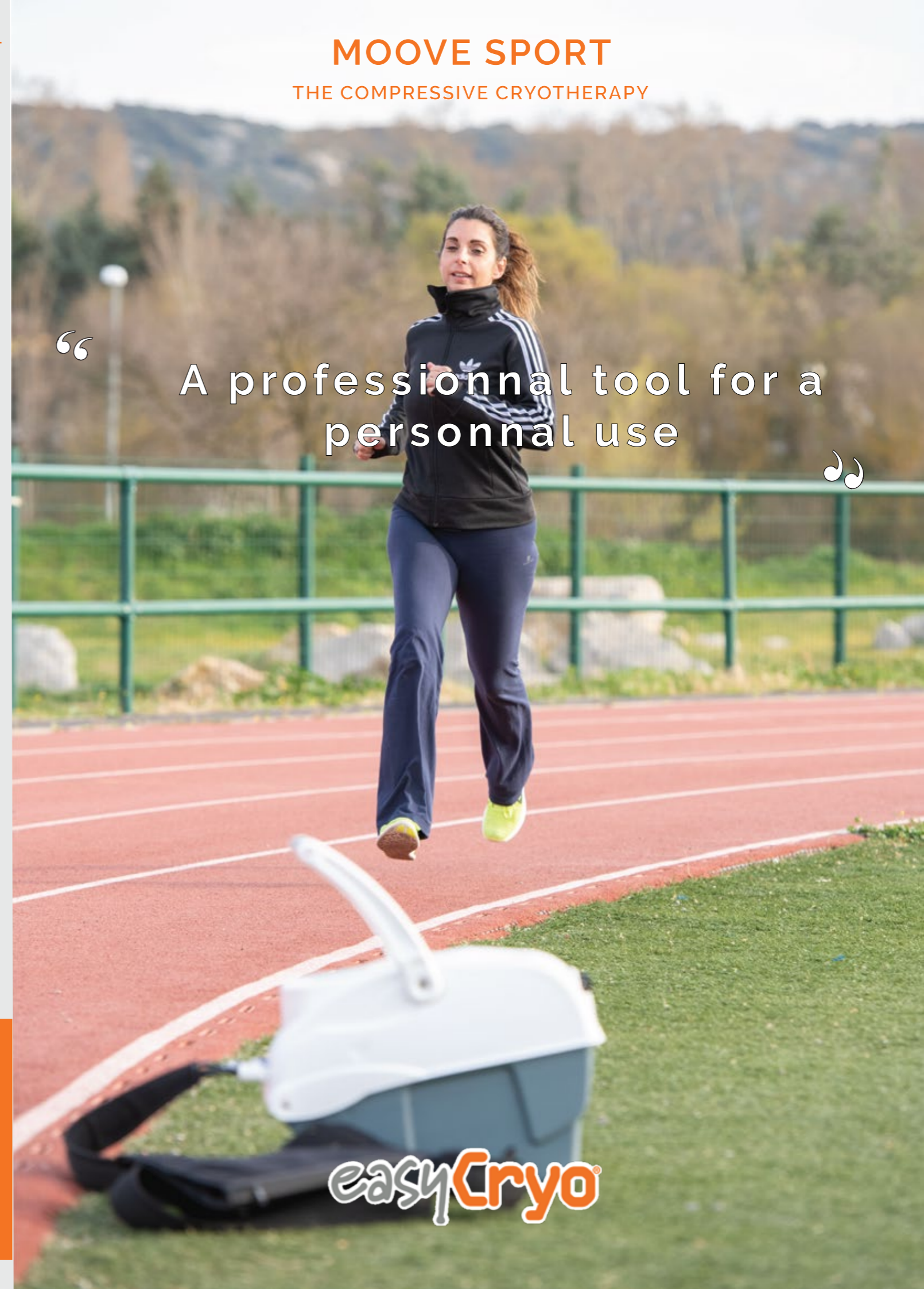

easyCryo®

Zone ECOPARC
27 Rue Mezière Christin
34690 FABRÈGUES - France



Tél. +33 (0)4 99 64 21 05
Fax. +33 (0)4 99 64 21 06
olivia.beldame@easycryo.fr



www.easycryo.fr

easyCryo®